

Single molecule force spectroscopy

#### AFM: lab on a tip

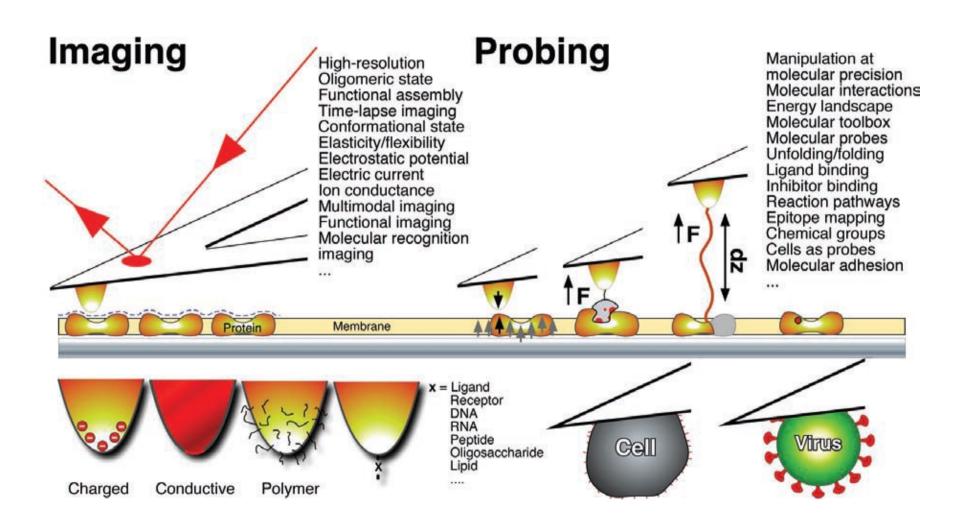


S. Scheuring, D. Muller, H. Stalhberg, H.-A. Engel, A. Engel, *Eur. Biophys. J.* **31**, 172 (2002)

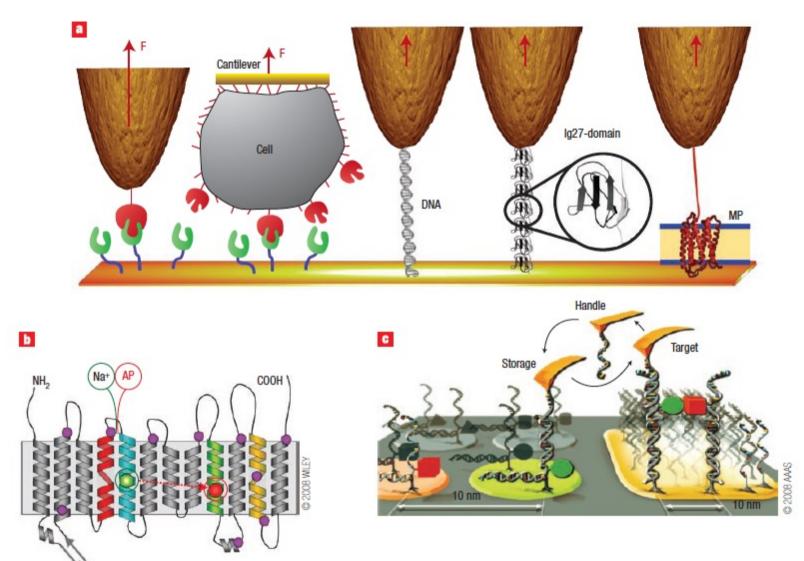


**Figure 3** Single molecule manipulation, control and design. **a**, Applying AFM to probe interaction forces (F) of single biomolecules. These examples measure ligand-receptor interactions in their isolated form (left) and embedded in their cellular environment (probe replaced by a biological cell); stretching of a DNA molecule; unfolding of Ig27-titin; and unfolding of a membrane protein (MP). **b**, SMFS can detect and locate interactions (circles) on the structure of membrane proteins (here proton/sodium antiporter NhaA from *E. coll*). A ligand (Na<sup>+</sup> ion) or inhibitor (AP, 2-aminoperimidine) binding to the ligand-binding site (green circle) establishes different interactions activating (green circle) or deactivating (green and red circles) the antiporter (composed of 12 transmembrane helices). **c**, Single molecules can be mechanically assembled by picking up from discrete storage sites with a DNA oligomer at the AFM probe and depositing them at a target site with nanoscopic precision. Reproduced with permission from refs 52, 53 and 63.

#### Single molecule force spectroscopy


Reversible Extension of a Single Dextran Polymer

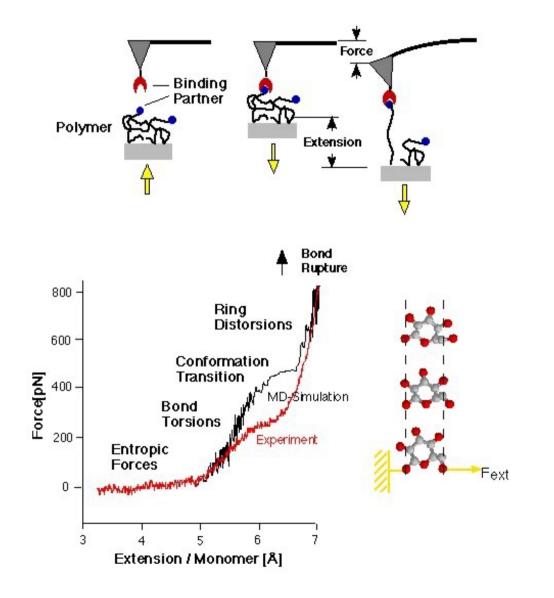

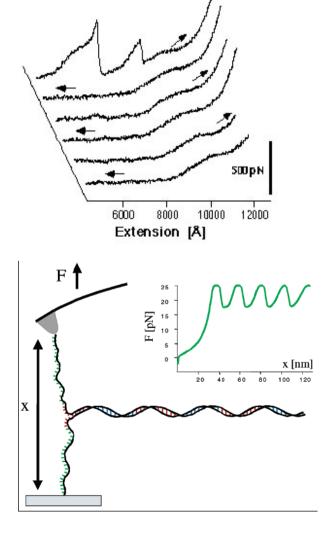


One common feature of many mechanical proteins is that they contain multiple, individually folded protein domains.

Two important examples are the **immunoglobulin (Ig)-type fold** and the **fibronectin-type fold** (the most common of which is fibronectin type 3 or FN-III). Both are so-called **b-sandwich structures:** these domains might unfold and refold as proteins execute mechanical functions.

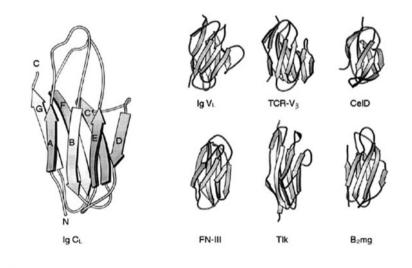

Force-induced extension of the **protein titin**, for example, which is responsible for the passive elasticity of muscle, can cause its constituent Ig and FN-III domains to unravel.


 $40 \times 25 \times 25 \text{ Å} (1 \text{ Å} = 10^{-10} \text{ m})$ 

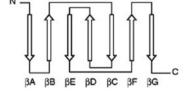


Structural diversity of Ig-type domains. The constant domain of antibodies (Ig  $C_L$ ) represents the common core present in all Ig-like structures, composed of seven antiparallel  $\beta$  strands labeled A–F. The two-dimensional topological diagram (bottom, left) shows the connected Greek key motif of this structure.

One common feature of many mechanical proteins is that they contain multiple, individually folded protein domains.

Two important examples are the **immunoglobulin (Ig)-type fold** and the **fibronectin-type fold** (the most common of which is fibronectin type 3 or FN-III). Both are so-called **b-sandwich structures:** these domains might unfold and refold as proteins execute mechanical functions.

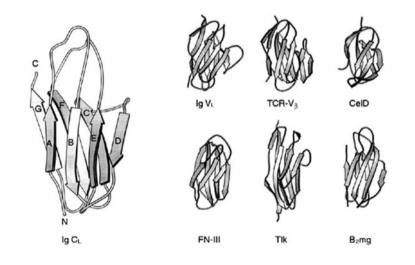


Unfolding and refolding of domains, as part of a particular mechanical function, could be **a mechanism by which tension is maintained as a protein is extended or relaxed**. Unfolding might also contribute to the function of fibronectin by exposing cryptic protein interaction sites that are important in ECM assembly.

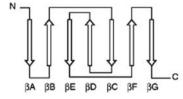

Structural diversity of Ig-type domains. The constant domain of antibodies (Ig C<sub>L</sub>) represents the common core present in all Ig-like structures, composed of seven antiparallel  $\beta$  strands labeled A–F. The two-dimensional topological diagram (bottom, left) shows the connected Greek key motif of this structure.


Single-molecule techniques analyze individual members in complex, heterogeneous populations, enabling to reveal the diversity of hidden molecular behaviors and rare events.

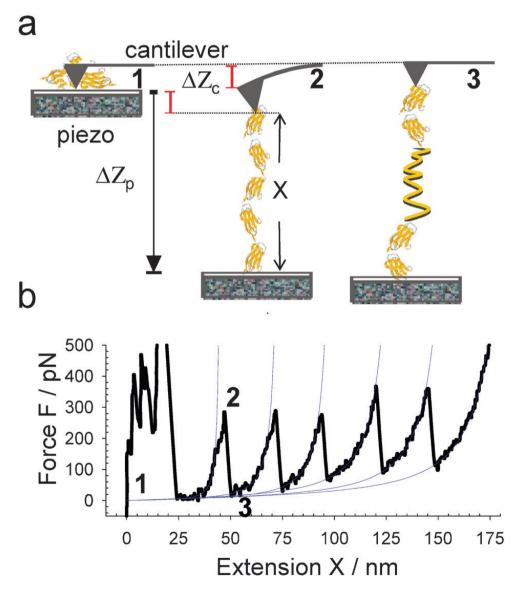

Single-molecule assays permit the direct measurement of the molecular elasticity of polysaccharides, proteins or nucleic acids, otherwise impossible, and the acceleration of certain molecular reactions, such as protein unfolding or ligand-receptor bond rupture.

$$X = \Delta Z_{\rm P} - \Delta Z_{\rm c}$$
$$F = k_{\rm c} \Delta Z_{\rm c}$$

P.E. Marszalek, Y. F. Dufrene Chem. Soc. Rev., 2012, 41, 3523–3534

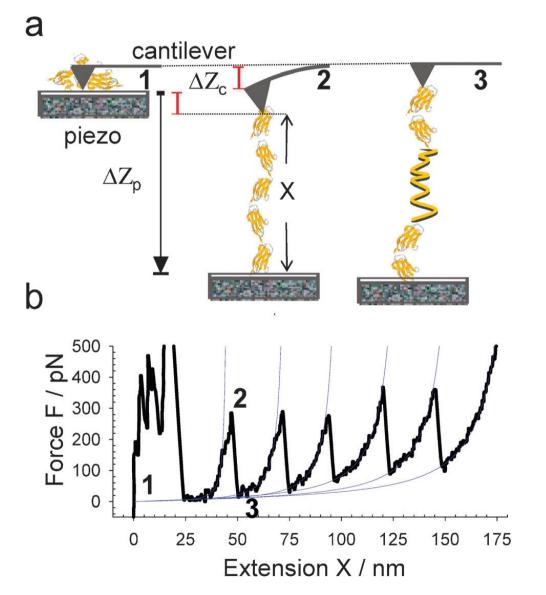

The resolution of force measurements depends on the dimensions of the AFM cantilever (the smaller the better) and on the range of measured frequencies

using small cantilevers in the measurement bandwidth of 1 kHz, the force noise that limits the resolution can be as small as 3–5 pN (resolution can be improved to below 1 pN)

 $\begin{aligned} X &= \Delta Z_{\rm P} - \Delta Z_{\rm c} \\ F &= k_{\rm c} \Delta Z_{\rm c} \end{aligned}$ 

#### How much should be the force to "break" a bond?

The energy of a molecular bond is approx 3 eV = 360 KJ / mole = 6 10(-19) J / bond

The length of a molecular bond is approx 0.3 nm or 3 Angstroem but to break a bond is sufficient to stretch it by less. Let us say 1 A

#### ENERGY = FORCE x DISTANCE


FORCE = ENERGY / DISTANCE = 6 10(-19) J / 0.1 10(-9) m

= 600 pN or less than 1 nN or 10 (-9) Newton

a cantilever  $\Delta Z_{c}$ 3 piezo  $\Delta Z_{p}$ Х b 500 Nd 400 400 400 2 Force | 200 100 0 0 25 75 100 125 150 175 50 Extension X / nm

Mechanical proteins are frequently composed of either individually folded globular domains that are connected to each other by short flexible linkers (e.g. titin in muscle) or tightly stacked short helical segments that form elongated single-domain structures (e.g. ankyrin in the RBC skeleton), or longer helical domains that form coiled-coils structures (e.g. spectrin in the RBC skeleton, or myosin in muscle).


P.E. Marszalek, Y. F. Dufrene Chem. Soc. Rev., 2012, 41, 3523–3534

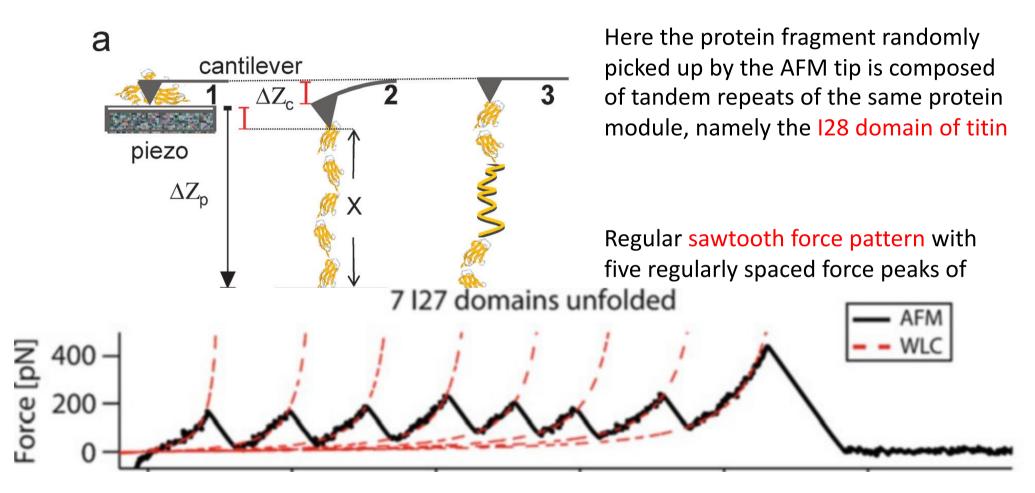

**Fig. 33.5** Typical force–extension plot of the unfolding of a polyprotein consisting of seven I27 domains from the titin protein (also called I91 domains). Each peak corresponds to an unfolding event of a single domain. The unfolding force for each domain is  $\sim$ 200 pN. The dashed red line indicates a family of worm-like chain fits with a contour length spacing of 28.5 nm between unfolding events


#### Polyprotein force spectroscopy

#### **Engineered poly I27 Titin**

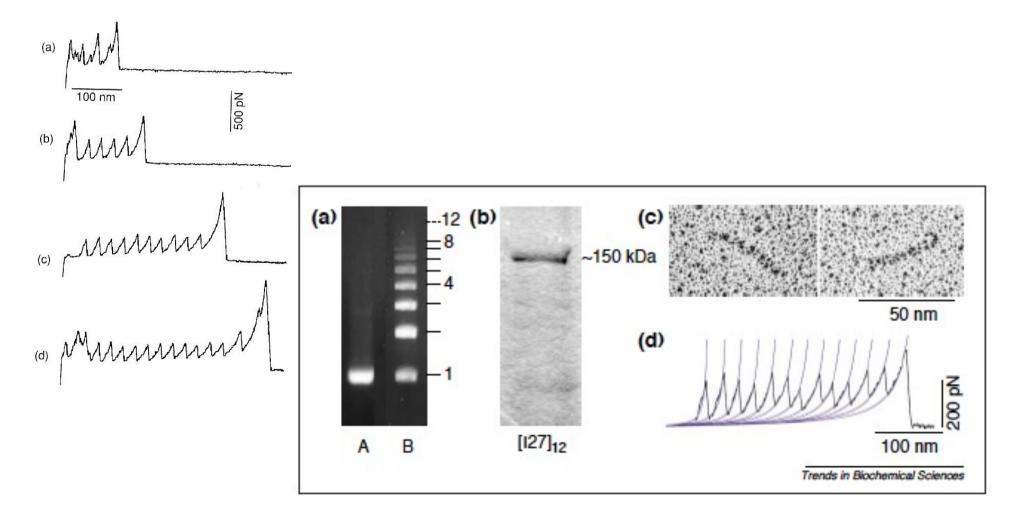


### Protein force spectroscopy

When a polymer is relaxed, it forms a coiled structure (maximizes the entropy of its segments). Extension of the polymer generates an opposing force (reduction in entropy) called **ENTROPIC ELASTICITY**: small extensions require little force, resistance to extension rises rapidly as the polymer approaches its full length

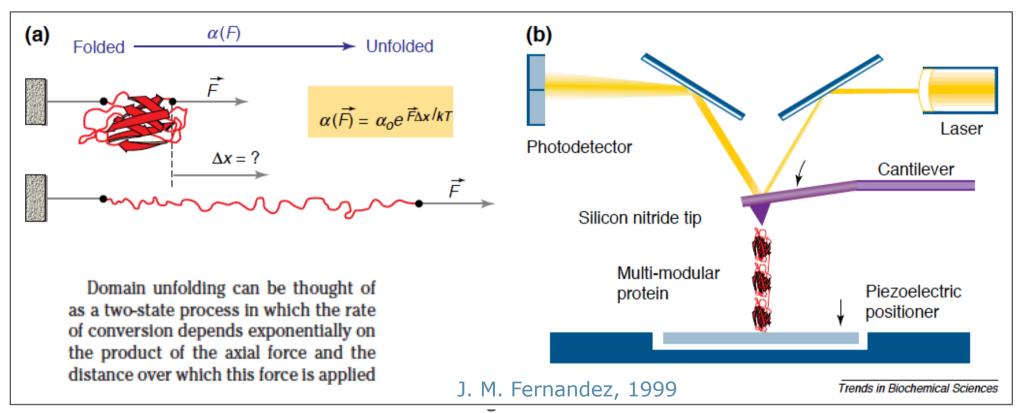

The unfolding of protein domains by an external force. (a) When axial stress is applied to a folded domain the protein will unravel. The inset shows an equation describing this transition, where *F* is the applied force,  $\Delta x$  is the distance over which the unfolding event occurs,  $\alpha_0$  is the rate constant in the absence of an applied force, *k* is Boltzmann's constant and *T* is the absolute temperature. Thus, the rate at which protein unfolding occurs increases exponentially with the applied force. This equation is similar to that describing the dissociation of non-covalent bonds placed under an external force<sup>38,39</sup>. (b) The force–extension mode of the atomic force microscope (AFM). When pressed against a layer of protein attached to a substrate, the silicon nitride tip can adsorb a single protein molecule. Extension of a molecule by retraction of the piezoelectric positioner results in deflection of the AFM cantilever. This deflection changes the angle of reflection of a laser beam striking the cantilever, which is measured as the change in output from a photodetector.


#### Polyprotein force spectroscopy

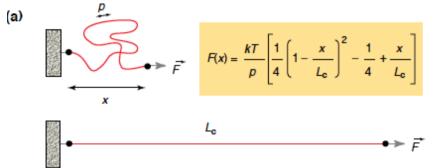


(b)


At small displacements, reduction in the number of conformations gives rise to entropic elasticity forces At large extensions, tension in the molecular backbone may lead to enthalpic elasticity effects (bond deformation, rupture of intramolecular hydrogen bonds and even conformational changes of the entire molecule)

The entropic elasticity of proteins and domain unfolding. (a) The entropic elasticity of proteins can be described by the WLC (worm-like chain) equation (inset), which expresses the relationship between force (F) and extension (x) of a protein using its persistence length (p) and its contour length ( $L_c$ ). k is Boltzmann's constant and T is the absolute temperature.

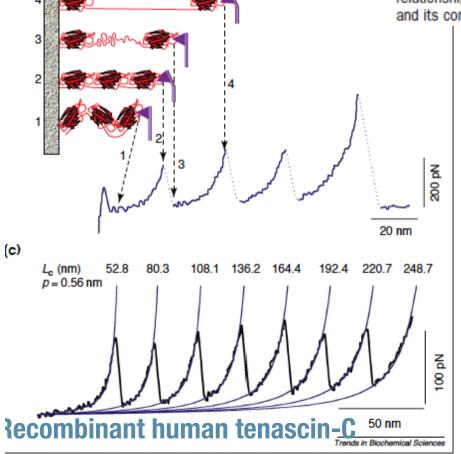

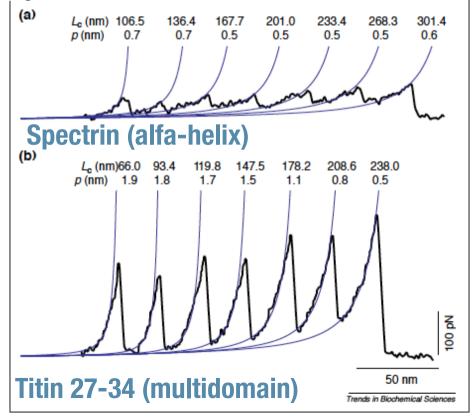

The behaviour of polymers under mechanical stress is described by the **wormlike chain (WLC) model of elasticity** 

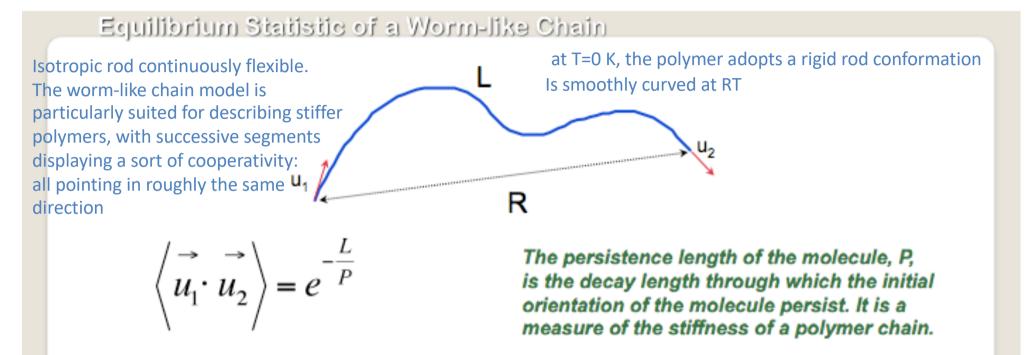

This model describes a polymer as a continuous string of a given total (or contour) length. Bending of the polymer at any point influences the angle of the polymer for a distance, referred to as the **persistence length**, that reflects the polymer flexibility. The smaller the persistence length, the greater the entropy of the polymer and the greater the resistance to extension.

The persistence length and the contour length comprise the adjustable parameters of the WLC model.

Courtesy of Prof. C. Bustamante

### Single molecule force spectroscopy

The statistical model used for describing the elastic behavior of the polymer is the WORM-LIKE CHAIN (WLC) model (proteins, DNA/RNA) and the FREELY-JOINTED CHAIN (FJC) model (polysaccharides)

#### WLC model

the polymer is an irregular curved filament, linear on the scale of the persistence length  $I_p$ , which represents the stiffness of the molecule. Molecules with low persistence length have a tendency to form coils. Extension is limited by the contour length  $L_c$  of the molecule, i.e. the length of the linearly extended molecule without stretching the molecular backbone.

In this model, the force **F versus** extension **x** is approximately given by:

$$F(x) = k_{\rm b}T/l_{\rm p}[0.25(1 - x/L_{\rm c})^{-2} + x/L_{\rm c} - 0.25]$$

## Single molecule force spectroscopy

The statistical model used for describing the elastic behavior of the polymer is the WORM-LIKE CHAIN (WLC) model (proteins, DNA/RNA) and the FREELY-JOINTED CHAIN (FJC) model (polysaccharides)

#### **FJC model**

the polymer is considered as a series of rigid, orientationally independent statistical (Kuhn) segments, connected through flexible joints. The segment length, or Kuhn length  $I_k$ , is a direct measure of the chain stiffness and is related to the contour length by  $Lc=nI_k$ .


In this model, the extension **x** versus force **F** is approximately given by:

$$x(F) = L_{\rm c}[\coth(Fl_{\rm k}/k_{\rm b}T) - k_{\rm b}T/Fl_{\rm k}]$$

And, allowing Kuhn segments can stretch and align under force with elasticity k<sub>s</sub>:

$$x(F) = L_{\rm c}[\coth(Fl_{\rm k}/k_{\rm b}T) - k_{\rm b}T/Fl_{\rm k}][1 + nF/k_{\rm s}L_{\rm c}]$$

#### Polyprotein force spectroscopy

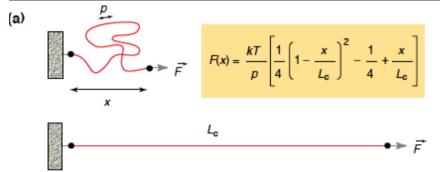



(b)

Interestingly, the measured unfolding forces generally increase with extension speed. Varying this parameter may provide an estimate of the unfolding off-rate and of the width of the unfolding energy barrier.


Reliable force data on a given molecule require recording several hundred force-curves using many independent tips and samples and proper model to fit data.

The entropic elasticity of proteins and domain unfolding. (a) The entropic elasticity of proteins can be described by the WLC (worm-like chain) equation (inset), which expresses the relationship between force (F) and extension (x) of a protein using its persistence length (p) and its contour length ( $L_c$ ). k is Boltzmann's constant and T is the absolute temperature.

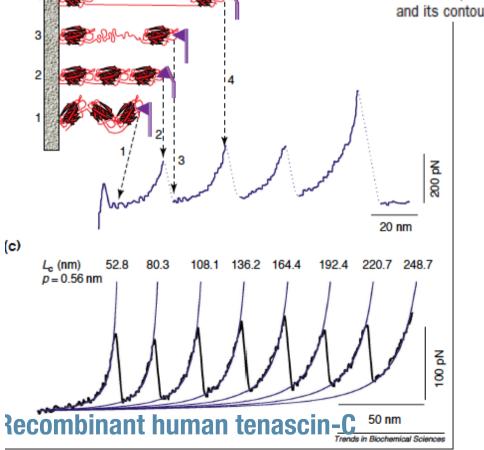

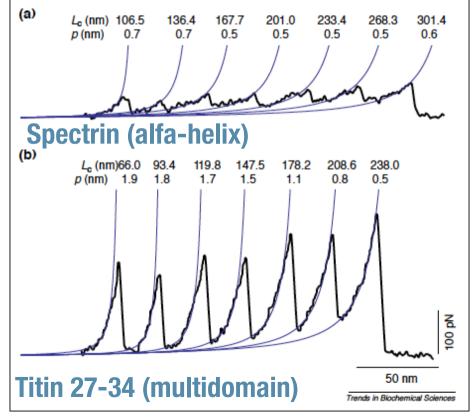


#### Unfolding-refolding

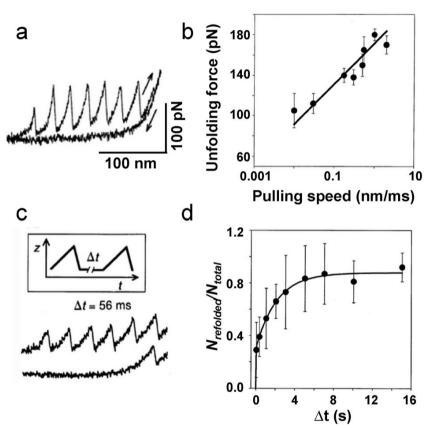

The data can be fitted using a Monte Carlo approach with a two-state model of the unfolding/refolding process, allowing the extraction of important parameters characterizing the unfolding reaction such as the distance to the unfolding transition state and the unfolding rate constant extrapolated to zero force

Fig. 3 Repeated unfolding and refolding cycles of tenascin. (a) The same TNfnAll fragment of tenascin was repeatedly stretched at various stretching speeds. (b) The average unfolding force vs. stretching speed. (c) A double-pulse experiment where the same molecule is stretched, relaxed for time  $\Delta t$  and stretched again. (d) The fraction of folded modules (counted from the second unfolding pulse) as a function of  $\Delta t$ .


# Polyprotein force spectroscopy

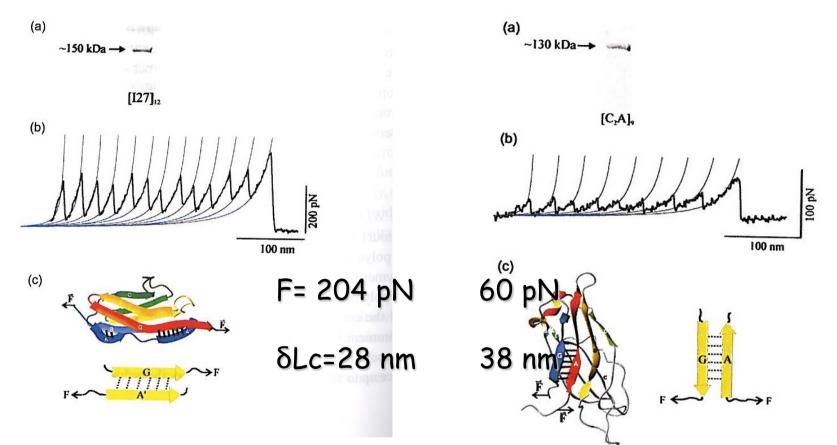

Two proteins with "all beta" structure and  $\beta$ -sandwich topology -titin has seven  $\beta$  strands which fold face-to-face through backbone Hbonds and hydrophobic core interactions, perpendicular to stretching direction

-C2  $\beta$  sandwich with 127 aminoacids arranged in 8 antiparallel strands.domains are in a zipper configuration


#### Herman E. Gaub (LMU Munich) : Unzipping DNA oligomers Nano Letters **3**, 4, 493 (2003)

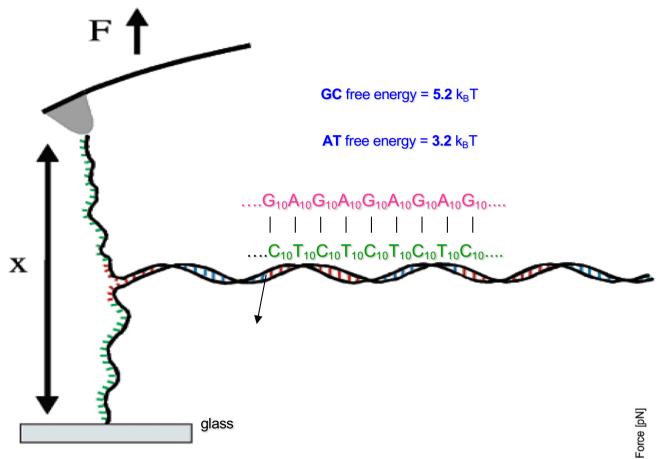



**Experimental Results** 

Complementary DNA oligonucleotides are chemically attached to an AFM tip and a glass slide. As the tip is brought in contact with the surface a ds-DNA forms, which can then be unzipped upon retraction the tip.

#### Herman E. Gaub (LMU Munich) : Unzipping DNA oligomers Nano Letters **3**, 4, 493 (2003)

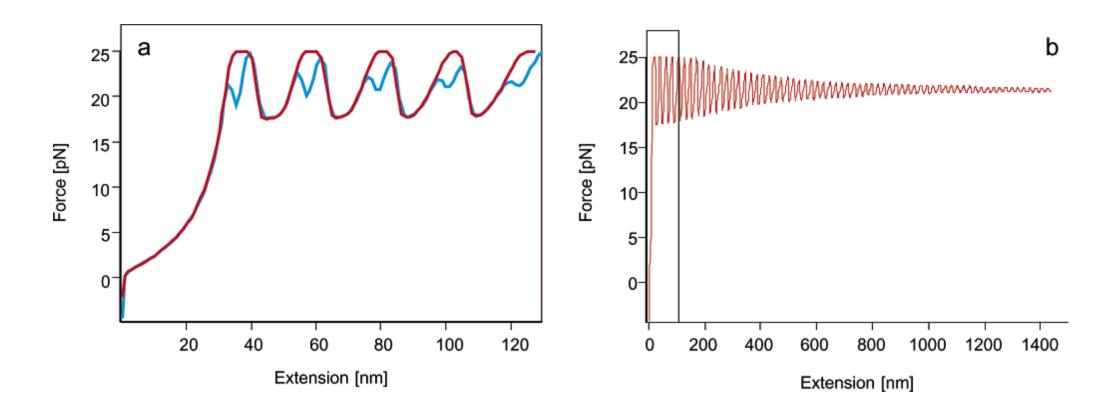

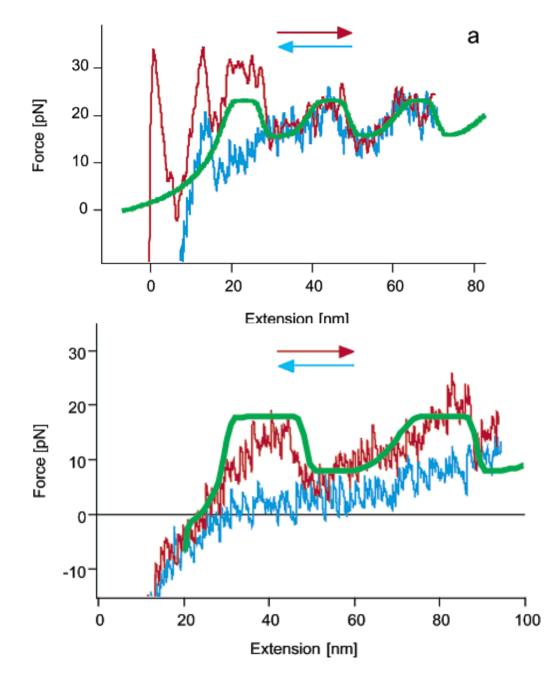

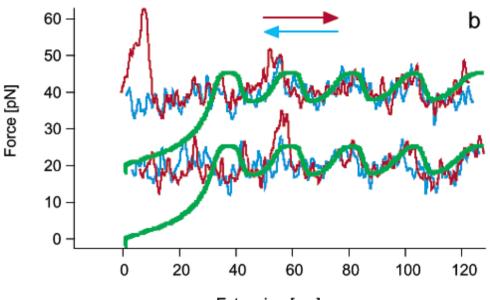

Extension [nm]

Figure 4. Force vs distance curve upon unzipping a dsDNA molecule with an alternating sequence of 20 pure AT and 20 pure GC base pairs ([dG20dA20dG20dA20dG20]([dC20dT20dC20dT20dC2 0]) and a simulated curve (green). The force upon unzipping the molecule varies because of the different stability of the AT and GC base pairs. The oscillation in force seen in the data corresponds well to the force vs distance profile calculated from the equilibrium thermodynamic model.